


**International Acceptance of the Nonradioactive LLNA: DA for Evaluating Allergic
Contact Dermatitis Hazards**

P Brown¹, J Matheson², A Jacobs¹, R Ward³, E Margosches³, E Salicru⁴, D Allen⁴, F
Stack⁴, W Stokes⁵

¹U.S. FDA, Silver Spring, MD, USA; ²U.S. CPSC, Bethesda, MD, USA; ³U.S. EPA,
Washington, DC, USA; ⁴ILS, Inc., RTP, NC, USA; ⁵NICEATM/NIEHS/NIH/DHHS, RTP,
NC, USA

Introduction


- The murine local lymph node assay (LLNA) is a test method for assessing the potential of substances to cause allergic contact dermatitis (ACD). ACD is an allergic skin reaction characterized by redness, swelling, and itching that can result from repeated contact with a sensitizing substance.
- In response to a nomination by the U.S. Consumer Product Safety Commission in 2007, the NICEATM evaluated the nonradioactive LLNA: DA (**Figure 1**) to assess the ACD hazard potential of substances.
- Daicel Chemical Industries, Ltd., developed the LLNA: DA (Yamashita et al. 2005; Idehara et al. 2008).
 - Measures ATP content in draining auricular lymph nodes as an estimate of cell number for the assessment of lymph node cell proliferation.
- ICCVAM published recommendations on the LLNA: DA in a test method evaluation report (available at: <http://iccvam.niehs.nih.gov/methods/immunotox/llna-DA/TMER.htm>).

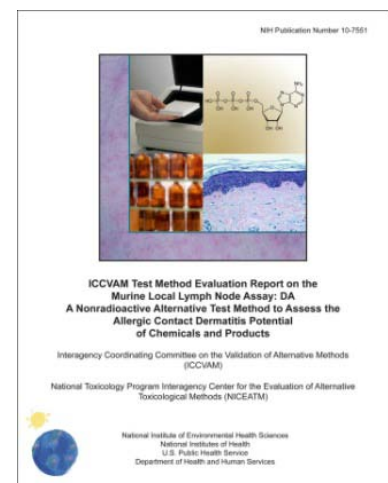

LLNA: DA Test Method Protocol

- The LLNA: DA protocol (Figure 1) is the same as the traditional LLNA protocol except for those procedures unique to the conduct of the LLNA: DA:
 - Pretreatment with 1% sodium lauryl sulfate prior to test substance application
 - One additional day of test substance application after 3 days of no treatment
 - Assessment of proliferation by measuring intracellular ATP levels within lymph node cells instead of ³H-thymidine incorporation
- The reduced LLNA: DA (rLLNA: DA) should be used routinely to determine the ACD hazard potential of chemicals and products in testing situations where dose-response information is not required or negative results are anticipated.
 - Like the reduced LLNA (Kimber et al. 2006; ESAC 2007; ICCVAM 2009), the rLLNA: DA protocol uses only the high dose and reduces animal use by up to 40%.
 - If existing information suggests a substance might have ACD hazard potential *and* dose-response information is needed, consider testing in the multidose LLNA: DA.

Figure 1 LLNA: DA Test Method Protocol

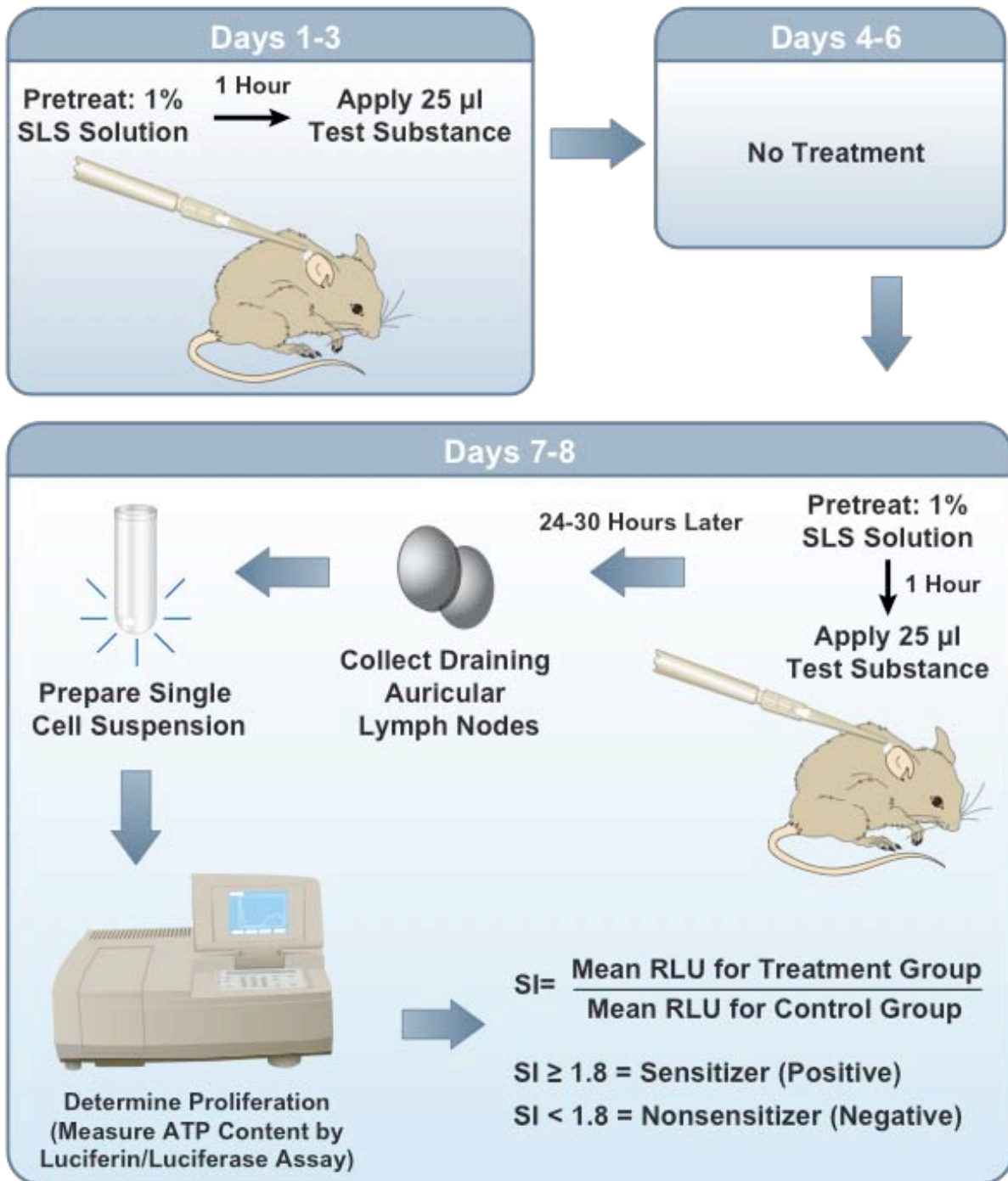

Abbreviations: RLU = relative luminescence units; SI = stimulation index; SLS = sodium lauryl sulfate

Current Validation Status of the LLNA: DA

Accuracy


- LLNA: DA database of 44 substances
 - Idehara et al. 2008
 - Idehara, unpublished data
 - Omori et al. 2008 (interlaboratory validation study)
- Results compared to traditional LLNA data
- Stimulation index (SI) ≥ 1.8 produced optimal results based on no false negatives
(Figure 2)
- The LLNA: DA correctly identified all 32 LLNA sensitizers (0% [0/32] false negatives) and 9/12 LLNA nonsensitizers
 - Accuracy = 93% (41/44)
 - False positive rate = 25% (3/12)
 - Chlorobenzene, hexane, and salicylic acid: all $1.8 < SI < 2.5$
 - False negative rate = 0% (0/32)

Reliability

- A concordance analysis of sensitizer (10/14) and nonsensitizer (4/14) outcomes was conducted across two phases of an interlaboratory validation study.
 - Concordance was observed for 80% (8/10) of the sensitizer outcomes.

- Two LLNA sensitizers, 3-aminophenol (1/3 SI < 1.8 and 2/3 SI ≥ 1.8) and nickel (II) sulfate hexahydrate (4/8 SI < 1.8 and 4/8 SI ≥ 1.8) produced discordant LLNA: DA test results.
- Concordance was observed for 75% (3/4) of the nonsensitizer outcomes.
 - The discordant LLNA nonsensitizer was isopropanol (91% concordance).

Figure 2 SI Decision Criteria Performance of the LLNA: DA Compared with the Traditional LLNA Using 44 Substances

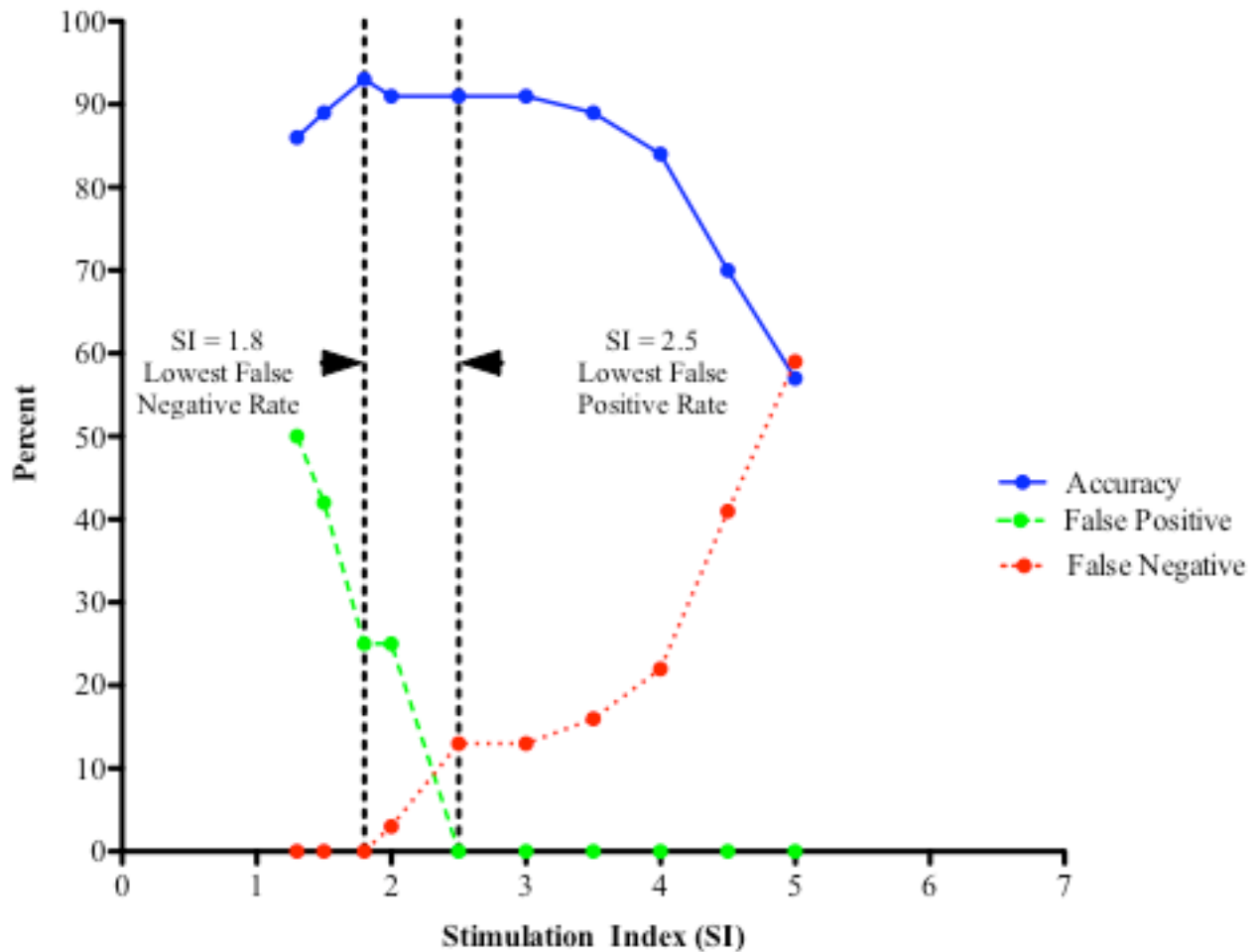


Compared to traditional LLNA results, the lines show the change in performance characteristics for the LLNA: DA with the SI used to identify skin sensitizers. This analysis used LLNA results for 32 sensitizers and 12 nonsensitizers. For 14 substances with multiple LLNA: DA test results, the most prevalent outcome was used.

LLNA: DA Test Method Usefulness and Limitations

- The LLNA: DA can be used to identify potential skin sensitizers or nonsensitizers.
 - Use $SI \geq 1.8$ to identify potential skin sensitizers.
- A slight potential for false positives with borderline weak positive responses ($1.8 < SI < 2.5$) exists.
 - Consider additional information such as dose-response relationship strength, statistical significance, evidence of systemic toxicity and/or excessive skin irritation together with SI values.
- The LLNA: DA might not be appropriate for testing certain classes of materials with properties that interfere with the assay. Consider if test substance might affect:
 - ATP levels (e.g., ATP inhibitors)
 - Accurate intracellular measurement of ATP levels (e.g., ATP-degrading enzymes or extracellular ATP in the lymph node)

LLNA Peer Review Panel Meetings

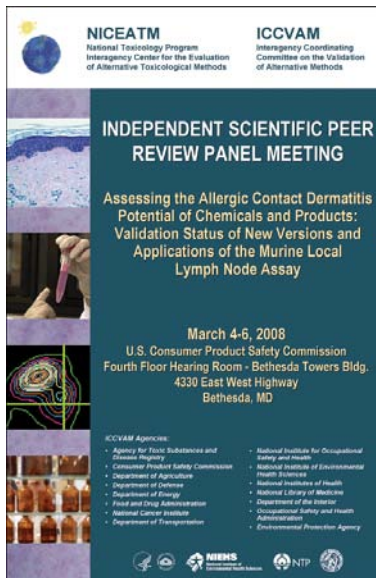

- Public meetings of an international independent scientific peer review panel were held at the Consumer Product Safety Commission in Bethesda, MD, on March 4-6, 2008, and at the National Institutes of Health in Bethesda, MD, on April 28-29, 2009.

Charge to the Peer Review Panel

- Review the draft Background Review Document (BRD) for errors and omissions
- Provide conclusions and recommendations on the current validation status of the LLNA: DA
- Does the information contained in the draft BRD support ICCVAM's draft test method recommendations?

Peer Review Panel Conclusions

- Concurred that the available data and test method performance supported the use of the LLNA: DA to identify substances as sensitizers and nonsensitizers, with certain limitations
- Concurred that, before additional animal testing is conducted, consideration should be given to the necessity for the substance to be tested for skin sensitization potential
- The complete LLNA Peer Review Panel Reports can be accessed at:
 - http://iccvam.niehs.nih.gov/docs/immunotox_docs/LLNAPRPRpt2008.pdf
 - http://iccvam.niehs.nih.gov/docs/immunotox_docs/LLNAPRPRpt2009.pdf

Independent Scientific Peer Review Panel

Michael Luster, Ph.D. (Panel Chair)

Senior Consultant to the National
Institute for Occupational Safety and
Health
Morgantown, WV

Nathalie Alépée, Ph.D.

L'Oréal Research and Development
Aulnay sous Bois, France

Anne Marie Api, Ph.D.

Research Institute for Fragrance
Materials
Woodcliff Lake, NJ

Nancy Flournoy, M.S., Ph.D.

University of Missouri–Columbia
Columbia, MO

Thomas Gebel, Ph.D.

Federal Institute for Occupational Safety
& Health
Dortmund, Germany

Sidney Green, Ph.D.

Howard University
Washington, D.C.

Kim Headrick, B.Admin., B.Sc.

Health Canada
Ottawa, Ontario, Canada

Dagmar Jírová, M.D., Ph.D.

National Institute of Public Health
Prague, Czech Republic

David Lovell, Ph.D.

University of Surrey
Guildford, Surrey, U.K.

Howard Maibach, M.D.

University of California–San Francisco
San Francisco, CA

James McDougal, Ph.D.

Wright State University
Dayton, OH

Michael Olson, Ph.D.

GlaxoSmithKline
Research Triangle Park, NC

Raymond Pieters, Ph.D.

Utrecht University
Utrecht, The Netherlands

Jean Regal, Ph.D.

University of Minnesota Medical School
Duluth, MN

Jonathan Richmond, MB ChB, FRCSEd

Home Office
London, U.K.

Peter Theran, V.M.D.

Consultant, Massachusetts Society for
the Prevention of Cruelty to Animals
Novato, CA

Stephen Ullrich, Ph.D.

M.D. Anderson Cancer Center
Houston, TX

Michael Woolhiser, Ph.D.

Dow Chemical
Midland, MI

Takahiko Yoshida, M.D., Ph.D.

Asahikawa Medical College
Hokkaido, Japan



International Acceptance of the LLNA: DA

- ICCVAM agreed with the OECD Expert Consultation Group that a single SI ≥ 1.8 to classify substances as skin sensitizers would avoid false negative and indeterminate results, which are not useful for regulatory purposes.
- OECD Test Guideline 442A Skin Sensitization: Local Lymph Node Assay: DA, which includes the SI ≥ 1.8 to classify substances as skin sensitizers, was adopted on July 22, 2010 (OECD 2010).
- OECD Test Guideline 442A can be accessed at <http://www.oecd-ilibrary.org/>
- International acceptance of the LLNA: DA is expected to result in broader use of LLNA tests.
 - Will further reduce and refine animal use for ACD hazard assessments on a global basis, while ensuring human safety
 - Will reduce costs and environmental hazards associated with the use of radioactive substances

ICCVAM Interagency Immunotoxicity Working Group

Consumer Product Safety Commission

Joanna Matheson, PhD (Working Group Co-chair)

Marilyn Wind, PhD (to July 2010)

Environmental Protection Agency

Office of Pesticide Programs

Jonathan Chen, PhD

John R. "Jack" Fowle III, PhD, DABT

Masih Hashim, DVM, PhD

Marianne Lewis

Deborah McCall

Timothy McMahon, PhD

John Redden

Jenny Tao, PhD

Office of Pollution Prevention and Toxics

Elizabeth Margosches, PhD

Ronald Ward, PhD

Office of Research and Development

Marsha Ward, PhD

Food and Drug Administration

Center for Devices and Radiological Health

Vasant G. Malshet, PhD, DABT

Jeffrey Toy, PhD

Center for Drug Evaluation and Research

Ruth Barratt, PhD, DVM

Paul Brown, PhD

Abigail Jacobs, PhD (Working Group Co-chair)

Jiaqin Yao, PhD

Center for Food Safety and Applied Nutrition

Donnie Lowther

Neil Wilcox, DVM, MPH (to April 2011)

Office of the Commissioner

Suzanne Fitzpatrick, PhD, DABT

National Institute of Environmental Health Sciences

Warren Casey, PhD, DABT

Dori Germolec, PhD

William Stokes, DVM, DACLAM

National Institute for Occupational Safety and Health

B. Jean Meade, DVM, PhD

Paul D. Siegel, PhD

National Library of Medicine

Pertti Hakkinen, PhD

European Centre for the Validation of Alternative Methods - Liaison

Silvia Casati, PhD

Alexandre Angers, PhD

Japanese Center for the Validation of Alternative Methods - Liaison

Hajime Kojima, PhD

References

ESAC. 2007. ESAC statement on the reduced local lymph node assay. ECVAM, Joint Research Centre, Institute for Health and Consumer Protection.

ICCVAM. 2009. ICCVAM Test Method Evaluation Report. The Reduced Murine Local Lymph Node Assay: An Alternative Test Method Using Fewer Animals to Assess the Allergic Contact Dermatitis Potential of Chemicals and Products. Available at <http://iccvam.niehs.nih.gov/methods/immunotox/LLNA-LD/TMER.htm>

Idehara K, et al. 2008. J Pharmacol Toxicol Meth 58:1-10.

Kimber I, et al. 2006. Contact Dermatitis 54:181-185.

OECD. 2010. Test No. 442A: Skin Sensitization: Local Lymph Node Assay: DA. In: OECD Guidelines for the Testing of Chemicals, Section 4: Health Effects, OECD Publishing. Available at http://www.oecd-ilibrary.org/environment/test-no-442a-skin-sensitization_9789264090972-en.

Omori T, et al. 2008. J Pharmacol Toxicol Meth 58:11-26.

Yamashita K, et al. 2005. AATEX 11:136-144.

Acknowledgments

The Intramural Research Program of the National Institute of Environmental Health Sciences (NIEHS) supported this poster. Technical support was provided by ILS, Inc., under NIEHS contract N01-ES-35504.

This poster reflects the views of the authors. The views expressed above have not been reviewed or approved by the U.S. Consumer Product Safety Commission or any other U.S. Federal agency, and do not necessarily represent the official positions of any U.S. Federal agency.

Since the poster was written as part of the official duties of the authors, it can be freely copied. NICEATM and ICCVAM gratefully acknowledge the following individuals and institutions for submitting data to NICEATM for the LLNA: DA evaluation.

Kenji Idehara, Ph.D.
Daicel Chemical Industries, Inc.
Hyogo, Japan

Takashi Omori, Ph.D.
Kyoto University School of Public Health
Kyoto, Japan